



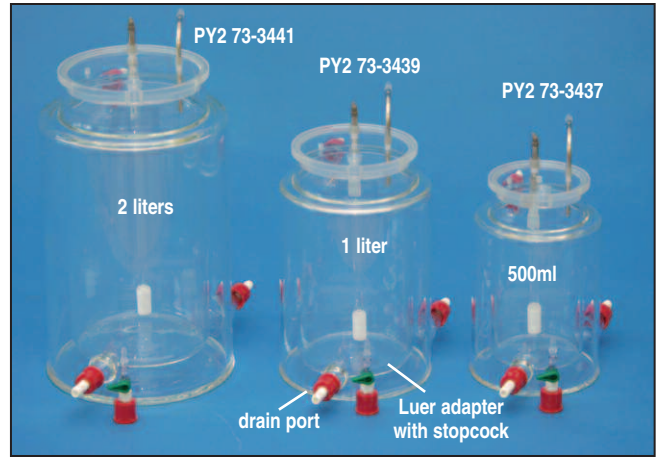

Glass Perfusate Reservoirs

accessories

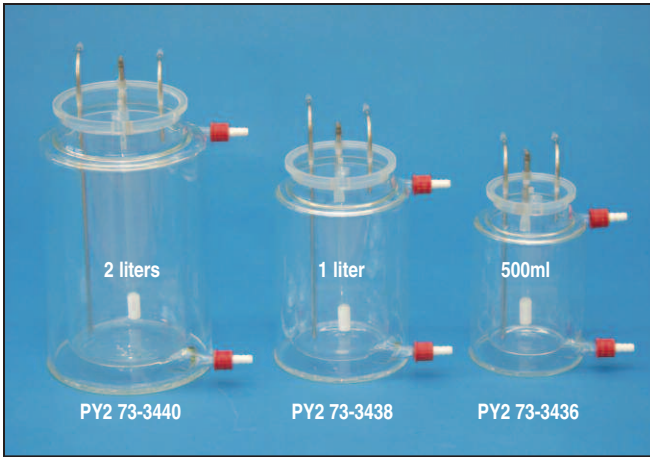

Glass Reservoirs

These Glass Reservoirs are for use with a pump. They are heated by a water thermocirculator. To connect to the thermocirculator 5mm ID tubing is necessary. For aeration of perfusate a large gas frit is included. The pump is connected to the longer stainless steel suction tube via a

LUER connector and tubing. A return flow can be connected to the short stainless steel tube with the included LUER Slip connector. Pictured tubing connections are included with each reservoir. Tubing sets must be purchased separately. See below.

Order #	Product
PY2 73-3440	2 Liter Reservoir Jacketed for Buffer Solution
PY2 73-3438	1 Liter Reservoir Jacketed for Buffer Solution
PY2 73-3436	0.5 Liter Reservoir Jacketed for Buffer Solution
PY2 73-0322	6 Liter Reservoir Jacketed for Buffer Solution with tube set for thermostating circuit and fluid line shutoff valves



Tubing Sets

Order #	Product
PY2 73-3455	Tube Set for Jacketed Buffer Reservoir
PY2 73-3456	Tube Set for Jacketed Buffer Reservoir with Fluid Line Shutoff Valves



PY2 73-T16428
3mm OD
Barbed Connector
(must be purchased separately)



PY2 73-3562
LUER Adapter
(included with the reservoirs above)

Reservoirs with Bottom Drain

These Glass Reservoirs have a bottom drain. They are for gravity supply systems or for use with a pump. The Reservoirs are heated by a water thermocirculator. For aeration of perfusate

they are supplied with a large gas frit. The outlet is the drain port at the bottom. It can be directly connected to a 5mm ID tubing or via a LUER adapter to the stopcock. A return flow can be connected to the short stainless steel tube via a LUER connector. Pictured tubing connections are included with each reservoir. Tubing sets must be purchased separately.

Order #	Product
PY2 73-3441	2 Liter Reservoir Jacketed for Buffer Solution with Bottom
PY2 73-3439	1 Liter Reservoir Jacketed for Buffer Solution with Bottom
PY2 73-3437	0.5 Liter Reservoir Jacketed for Buffer Solution with Bottom
PY2 73-3566	Repl. Frit for 2.0 Liter Glass Reservoirs
PY2 73-3565	Repl. Frit for 1.0 Liter Glass Reservoirs
PY2 73-3564	Repl. Frit for 0.5 Liter Glass
PY2 73-3562	Fluid Outlet with Stopcock for Jacket Reservoir with Bottom Drain



Connect With Us

