


Staining Machine MSM

The most universal unit easy to operate

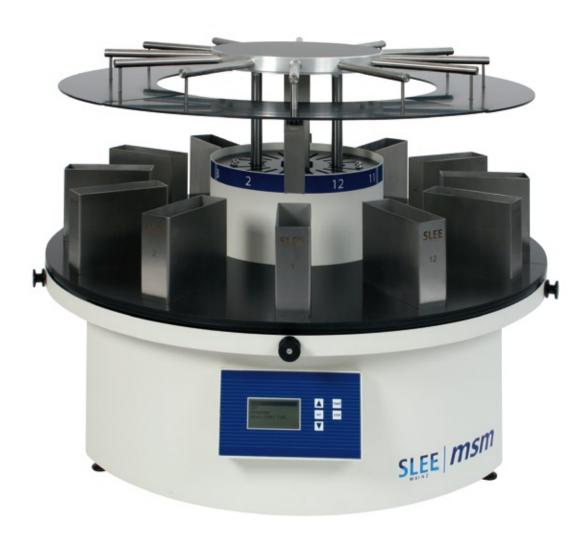

Staining Machine MSM

The most universal unit easy to operate

Universal unit with reduced staining time for better stain reproducibility. Usable for histology and cytology.

- Up to 1150 slides per hour
- 10 staining programs, free programmable
- ► Display and operating panel with programming keys
- ▶ LCD display with 4 lines and only 5 operating keys, free programmable
- ▶ Possibility to store up to 10 programs with individual rest and dripping time for all baths
- ► Display of active baths and remaining time
- Possibility of a start delay of up to 99 hours
- ► All exterior parts are powder coated in RAL 9002
- ► Cuvette locking made of black PVC
- 12 stainless steel cuvettes, 150 x 50 x 120 mm, 400 ml content
- 1 transport basket is included, its capacity is 40 standard size micro slides (20 or 50 optional), standard size 76 x 26 mm
- Rear connection possibility for hose (100 mm dia.)
- Electronic control with multi power supply 100 240 V, 50/60 Hz (worldwide operation without power adjustment)
- In case of power failure the basket is moved to a programmable position
- Stainless steel baskets for receiving the micro slides. They will be filled outside and then closed with a snap lock. Filled baskets can now be loaded to the receiving
- ▶ Lifting unit for lifting and lowering loaded baskets and closing lids available for each cuvette
- Rotation unit moves all cuvettes in programmed sequence
- Optional:
 - 24 cuvettes
 - Perspex hood for covering moveable equipment parts during operation and as suction hood for stopping bad odours, (side opening for container loading)
 - Exhaust fan
 - Exhaust fan carbon filter
 - Up to 7 transport baskets with continuous staining

Weight and Dimensions Dimensions (in mm) Dimensions packed (in mm)	Depth x Width x Height 870 x 870 x 700 1030 x 980 x 800
Weight net Weight packed	80 kg 140 kg
Electrical parameters	ArtNo. 12000200 230 V / 50Hz
	ArtNo. 12000201 230 V / 60Hz ArtNo. 12000202 115 V / 60Hz
	ArtNo. 12000203 100 V / 50Hz ArtNo. 12000204 100 V / 60Hz

Please ask for other units for Histopathology and consumables as well as laminar flow equipment.





MSM during staining



Cuvettes with transport basket



Transport basket for slides





Connect With Us











